
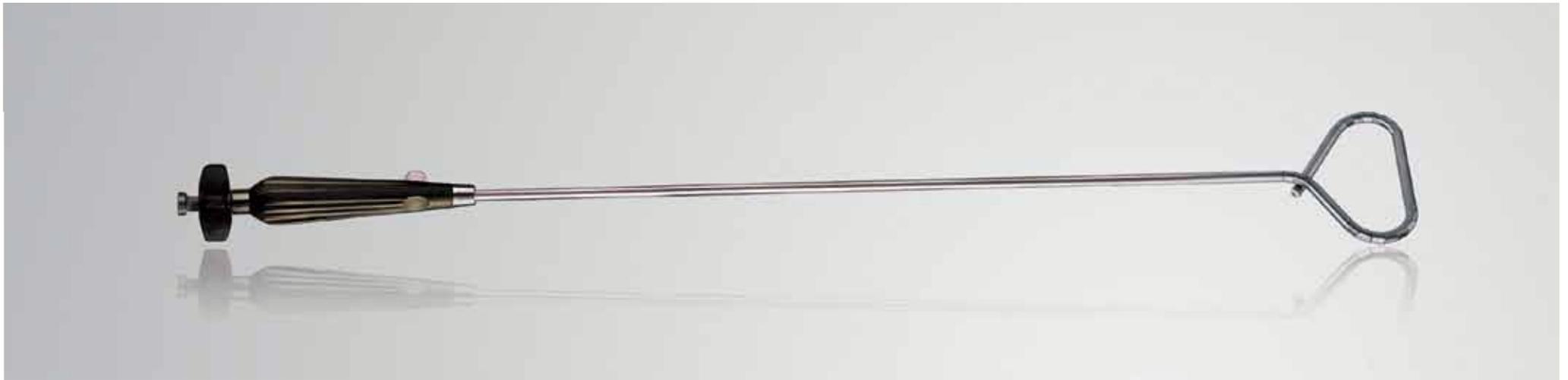



Diamond-Flex® laparoscopic articulating retractors

V. Mueller® and Snowden-Pencer® instrumentation

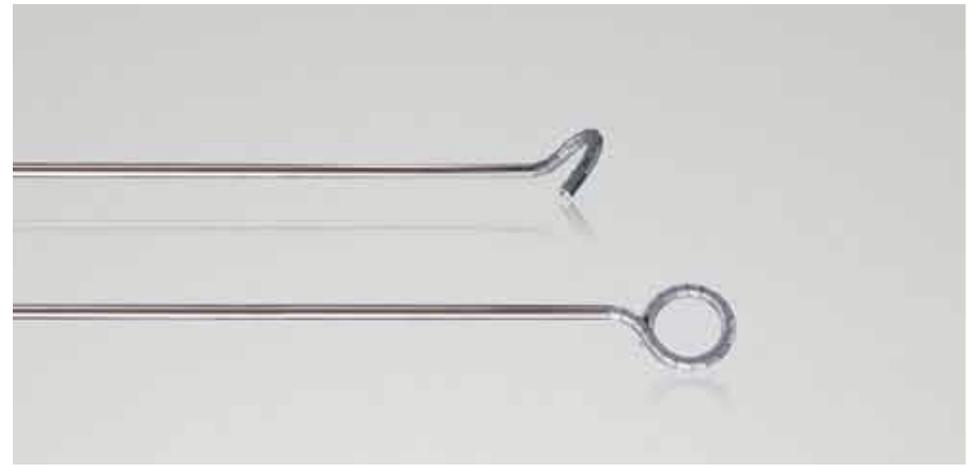

Diamond-Flex® technology allows instrumentation to articulate in the peritoneal space for retraction of organs and structures.


This exclusive line of reusable instruments helps meet the need of unique surgical challenges by enabling the surgeon to insert the instrument straight through a 5.5mm cannula and then activate one of many predetermined shapes. To activate, articulate the tip of the instrument into shape by turning the actuating knob on the handle until the complete shape is formed.

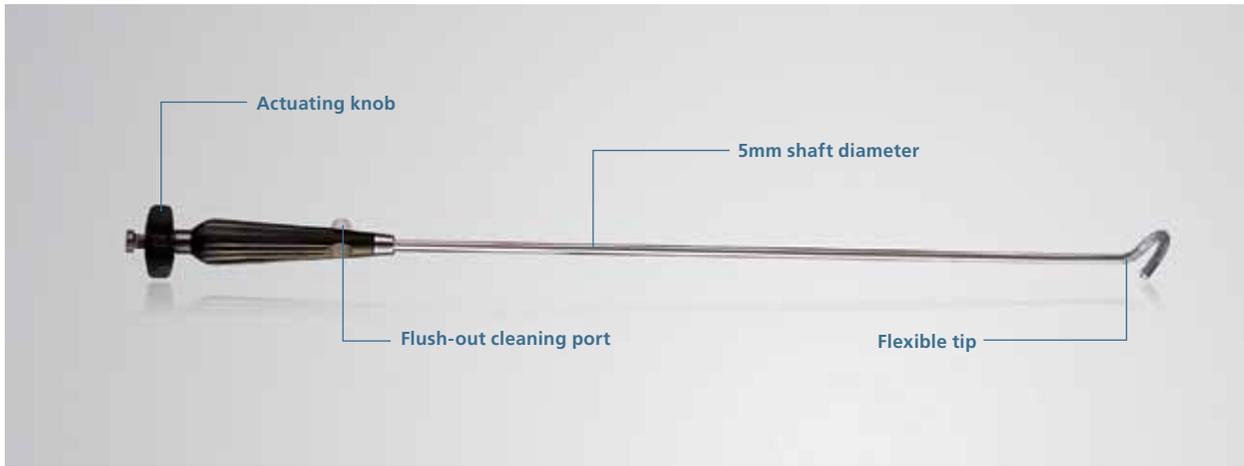

Catalog No.	Description	Shaft length
89-6101	Angled retractor	36cm
89-6104	Hooked retractor, curved right	36cm
89-6105	Hooked retractor, curved left	36cm
89-6106	Angled retractor , curved right	36cm
89-6107	Angled retractor, curved left	36cm
89-6109	60mm triangular retractor	36cm
89-6110	80mm triangular retractor	36cm
89-6210	80mm triangular retractor	45cm
89-6216	120mm triangular retractor	45cm
89-6111	60mm angled triangular retractor	36cm
89-6212	80mm angled triangular retractor	45cm
89-6112	80mm angled triangular retractor	36cm
89-6113	40mm circular retractor	36cm
89-6114	40mm angled circular retractor	36cm
89-6119	Protective sterilization sleeve	N/A
88-9861	Flush port plug	N/A

All shaft diameters are 5mm

Diamond-Flex® retractors have a flush port system to aid in cleaning and sterilization. In addition, the Diamond-Flex® retractor is accompanied by a metal sterilization sleeve used to keep the distal end straight during sterilization and cleaning.

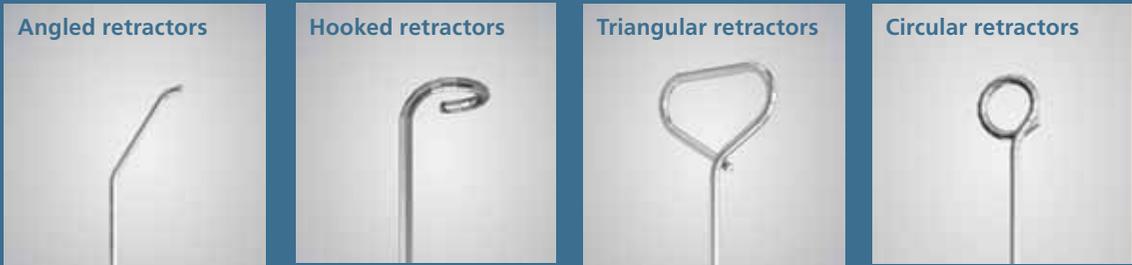

Snowden-Pencer® warranty

Diamond-Flex® retractors feature a one-year service warranty, which covers the most common repairs laparoscopic instruments require, and a life-time warranty against defects in materials and workmanship. Any product that proves to be defective in materials and workmanship will be repaired or replaced at our discretion. This guarantee (warranty) does not cover damage caused by abuse or misuse and is voided if anyone other than Snowden-Pencer® personnel repair or modify the instrument.

To place an order or for more information, please call your V. Mueller® and Snowden-Pencer® sales representative, call 800.323.9088 or visit carefusion.com/snowdenpencer.

© 2010 CareFusion Corporation or one of its subsidiaries. All rights reserved. Diamond-Flex, Snowden-Pencer and V. Mueller are registered trademarks of CareFusion Corporation or one of its subsidiaries. VM1561 (0410/PDF)